


## HemaPrep – Blood smearing instrument

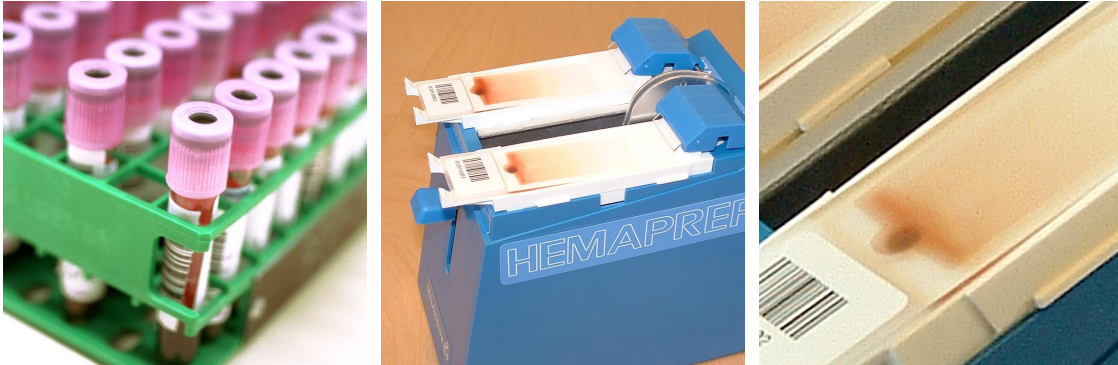

HemaPrep® is a blood smearing instrument that standardizes the making of smears. The instrument is small, portable and needs no power source. By adjusting the speed of the spreader, the operator is able to control smear length and thickness. HemaPrep is designed to prepare two slides at the time, from the same or different blood samples. HemaPrep has been available on the American market for over 15 years.

### Smear quality

A high quality smear is a condition for a reliable analysis. It takes a lot of training to make a high quality blood smear i.e. a smear of the right length and thickness. With HemaPrep the laboratory can standardize the smear quality irrespective of the staff experience.

### Easy to use

HemaPrep is very easy to use;

- Place the slide in position
- Apply 5-6µl of blood on the target area
- Depress the actuator bar carefully but firmly and release it
- Done

### Technical facts

*Recommended slide size:* 26x76mm.

*Instrument size:* approximately 90x135x155 mm (width x height x depth)

*Order number:* XU-10012